



## HEINE HL 1200® Examination Light



- ▷ Two flexible sections
- ▷ Ultra-slim head
- ▷ 90,000 Lux at 30 cm



**100% more light with new Bulb technology.** The HL 1200 is a flexible halogen examination light with a 20W IRC reflector lamp with innovative bulb technology. More than twice the brightness in comparison with conventional 20-Watt lamps for bright, white concentrated light. HL 1200 is a unique design incorporating a number of important features which most other lamps cannot match:

- ▷ **Ultra-slim head.** More freedom of movement and an unobstructed view for the examiner.
- ▷ **Built-in temperature filter.** Avoids patient discomfort.
- ▷ **Swivels through complete 180°.** Flexible and easy to position.
- ▷ **Exceeds all relevant quality standards.** Guaranteed protection for the patient even if e.g. a bulb fails.
- ▷ **“Stay-put” construction of flexible section.** Adjust the light to virtually any position and it “stays-put”, will not sag or droop.
- ▷ **Special non-tip wheeled stand.** Resists tipping to improve safety in clinic and practice.
- ▷ **Long-life 12V/20W IRC reflector lamp.** 4000 hours bulb life.
- ▷ **90,000 Lux light output at 30 cm.** Bright, white light.

### HL 1200 examination light, black

12V/20W

with wheeled stand (A), metal base	J-005.27.068
with lightweight wheeled stand (B), polycarbonate base	J-005.27.078
with additional multi-purpose clamp (C)	J-005.27.056
with wall-mounting bracket, can also be fitted to opt. wheeled stand (D)	J-005.27.052

### HL 1200 examination light, white

with wheeled stand (A), metal base	J-005.27.067
with lightweight wheeled stand (B), polycarbonate base	J-005.27.077
with additional multi-purpose clamp (C)	J-005.27.058
with wall-mounting bracket, can also be fitted to opt. wheeled stand (D)	J-005.27.054

Spare bulb 12V/20W IRC reflector	J-005.27.075
----------------------------------	--------------

Wheeled stand (metal base) with mount for HL 1200 light	Y-096.50.004
Wheeled stand (polycarbonate base) with mount for HL 1200 light	Y-096.50.009
Mount for HL 1200 light for wheeled stand	J-000.27.069



230V is the standard mains voltage. If a different voltage is required, e.g. 120V, please specify when ordering.

## HEINE HL 5000® Examination Light

### ► NEW DESIGN !



- ▷ Two flexible sections
- ▷ Exceptionally-high light output
- ▷ 140,000 Lux at 30 cm



With an exceptionally-high light output of 140,000 lux at a distance of 30 cm the HL 5000 examination light is ideal for examining and recognizing small details. The flexible gooseneck construction provides a working radius of 120 cm. Three versions are available: wheeled stand, wall mounting, multi-purpose clamp.

- ▷ **Long-lasting 12V/50W halogen bulb gives 140,000 Lux light output at 30 cm.** Very bright, white light.
- ▷ **Extremely compact illumination-head.** No obstruction of view, ease of handling.
- ▷ **3 different mounting options: wall, table or wheeled stand.** Versatile.
- ▷ **Halogen bulb lasts at least 4000 hours.** Economical.
- ▷ **Additional UV filter.** Increased safety during examination.
- ▷ **Solid construction with high-specification 120 cm radius gooseneck design.** Super-wide radius of operation, easy to position, light stays put in any selected position.
- ▷ **Additional IR-Filter reduces heat.** More comfort for the patient and examiner.

#### HL 5000 examination light

12V/50W

on non-tip wheeled stand (A), metal base

J-005.27.106

on lightweight wheeled stand L (B), polycarbonat base

J-005.27.108

with (C) wall mount

J-005.27.102

with (D) clamp mount

J-005.27.104

Spare bulb 12V/50W halogen

J-005.27.111



230V is the standard mains voltage. If a different voltage is required, e.g. 120V, please specify when ordering.



## HEINE HKL® Universal Fiber Optic Examination Light




The ideal light for any doctor's practice. Provides light similar to daylight. Can be powered by any HEINE F.O. projector.

- ✦ **Adjustable Diaphragm Aperture.** Wide range of spot sizes. Bright, even light. No temperature increase.
- ✦ **Colour temperature 3800K.** "Daylight" illumination.
- ✦ **Long, flexible arm "stays put" where required.** Very flexible to position with long-reach arm.
- ✦ **Slimline head.** Unobstructed view.
- ✦ **Universal-mounting.** Suits the practitioner's needs.
- ✦ **Can be retro-mounted to existing projectors.** Makes an examination light out of any projector.
- ✦ **30,000 Lux light output at 20 cm with projector HK 7000.**

Mounting options


- on wheeled stand
- with clamp for table-top mounting (maximum 40mm thick)
- with clamp for vertical mounting to a tube up to 40mm dia.
- with clamp for mounting to square tube up to 40mm dia.
- wall-mounting for HK 7000 projector.

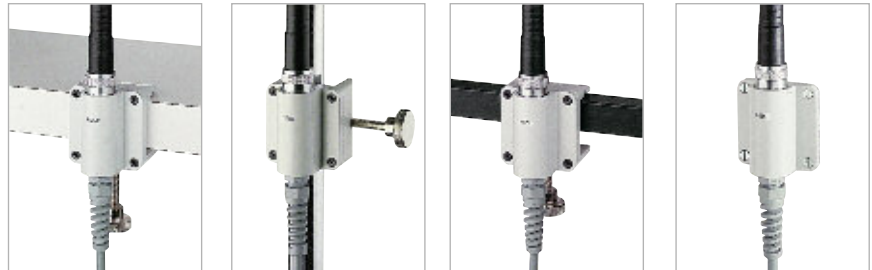

HKL Light	F.O.
with clamp	J-003.27.001
with wall-bracket	J-003.27.002
on wheeled stand	J-003.27.004
Clamp black, only	J-000.27.010
Clamp white, only	J-000.27.012
Base plate, only	J-000.27.011
Wheeled stand without base-plate	Y-096.50.001
Wheeled stand with base-plate	Y-096.50.002

Please order the projector separately

Light spot sizes

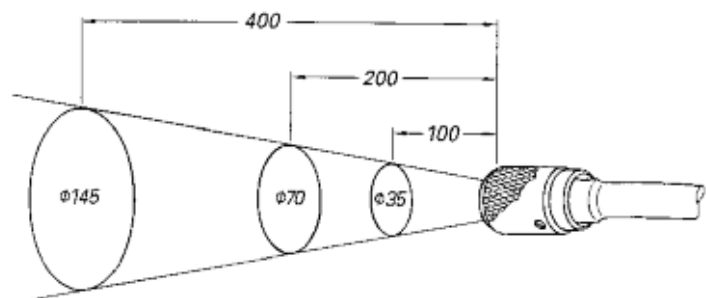

230V is the standard mains voltage. If a different voltage is required, e.g. 120V, please specify when ordering.