



RETINAL CAMERA

TRC-50DX

NEW



The Gold Standard in Retinal Imaging

•e IA), improves on the unsurpassed quality of Topcon retinal cameras, incorporating new functions that enhance their versatility and operational ease. The camera can be used with a variety of photo devices, from 35mm camera backs to super-high resolution digital cameras. A newly developed control panel, equipped with a touch-screen, allows for a quick and easy change of settings while the Small Aperture function facilitates focusing and enhances sharpness.

Further enhancing its adaptability, the TRC-50DX(Type IA) can perform Auto Fluorescence studies, expanding its diagnostics aptitudes and application availability.

An adaptable filter mechanism lets the user add or exchange filters of his or her choice for different procedures or research projects.


The TRC-50DX : Leading the way in retinal cameras.

RETINAL CAMERA **TRC-50DX / 50DX** (Type IA)

Features

Auto Fluorescence

This new diagnostic technique allows for a more in-depth study of the Retinal Pigment Epithelium (RPE) without the need of an injection of a contrast media, enhancing patient comfort and cooperation. This feature is available as the default component in TRC-50DX type IA.

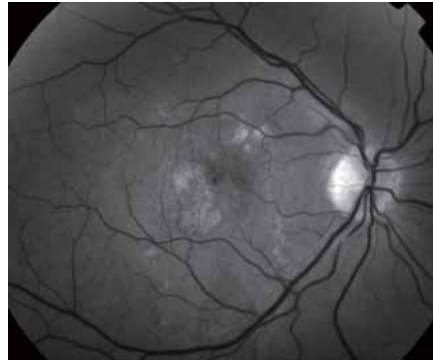


Red Free



Auto Fluorescence

Split Image Focusing

This convenient feature minimizes out-of-focus images and allows even less experienced operators to obtain clear, sharp images with minimal effort.

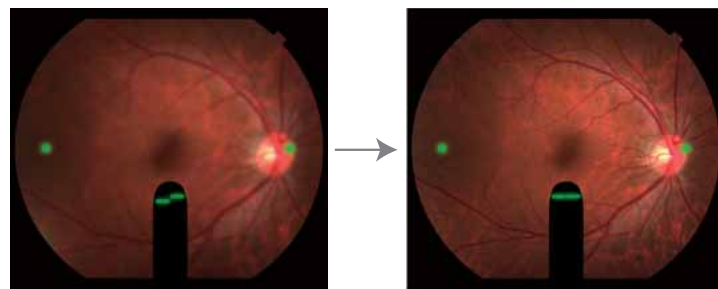

Compact Design

Minimized design helps to check patient's reaction and to open eye beyond the unit.

Flash Intensity Control

Conveniently illuminated for use in subdued light, the flash intensity control has 21 levels and can be used in full or half steps. Sensitive membrane buttons permit total control on the image brightness.



External Fixation Target

The newly designed fixation target guides the patient's eye to any desired position with a brilliant LED that shines in red or green light, changing automatically according to the study mode.



Touch Panel Control

Designed with user convenience in mind, the TRC-50DX features a Touch Panel control that allows the operator to change the camera settings with a single touch. All adjustments are at the photographer's fingertips.

Small pupil mode

The Small pupil mode permits imaging through pupils as small as 4.5mm diameter.

Aperture adjustment

The "Small" setting of the aperture produces a sharper, easier to focus image with enhanced depth of field. It offers remarkable effect in peripheral photography. Furthermore, its enhanced range of focus helps when using high resolution camera which requires severe focus aligning.




Off-focus





SMALL


peripheral

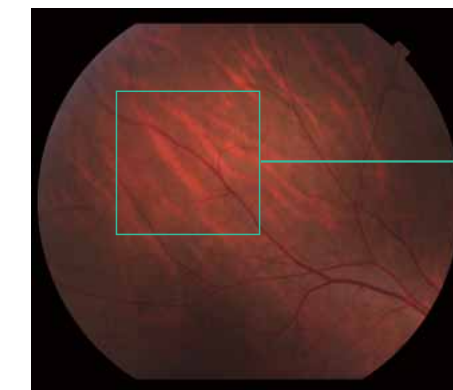

SMALL

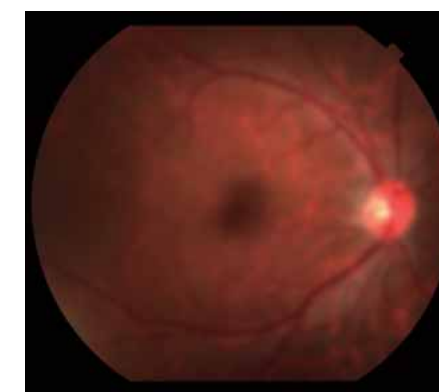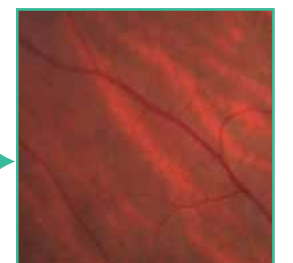

LARGE

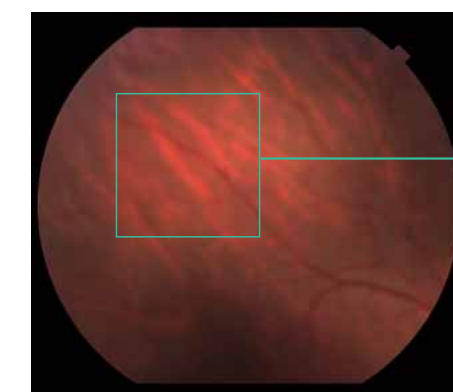

LARGE



Optional Accessories

RETINAL CAMERA TRC-50DX & TRC-50DX (Type IA)



Relay Lens

OR-2 (1X Relay Lens)

The OR-2 enables conventional color photography and fluorescein angiography on film.



Double Mount TV Relay Lens

Various double relay lenses meet the requirement of attaching CCD cameras to the top part.



Triple Mount TV Relay Lens

Triple mount TV relay lenses enable integration of 3 CCD cameras to the top part.







Single Mount TV Relay Lens

Single mount TV relay lenses make up simple and economical systems for color photography and fluorescein angiography.




IMAGeNet™ Digital Imaging System

Using comprehensive software for capturing, processing, archiving and consulting, IMAGeNet Professional expands the abilities of the TRC-50DX into a digital imaging analysis system for clinical applications. IMAGeNet Professional detects the filter settings and CCD camera in use from the TRC-50DX and sets the capture screen to the corresponding procedure automatically. A variety of options conveniently priced fit the needs of the single practitioner or the multi-discipline health institution with powerful diagnostic tools. The smart electronics of the TRC-50DX and the convenience of USB connections provide a seamless interface with any of the IMAGeNet Professional programs.

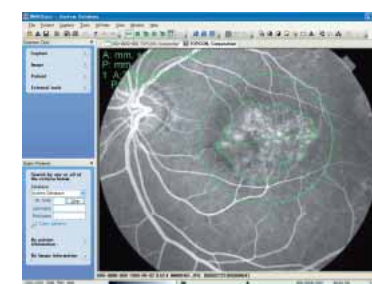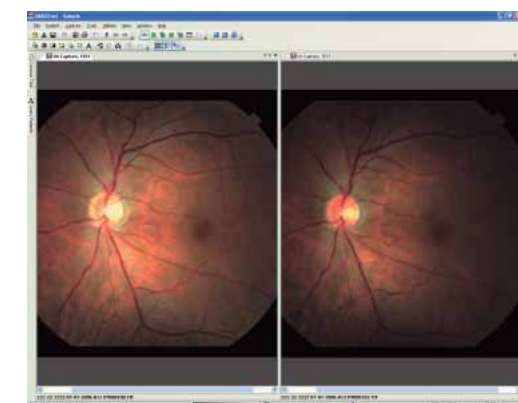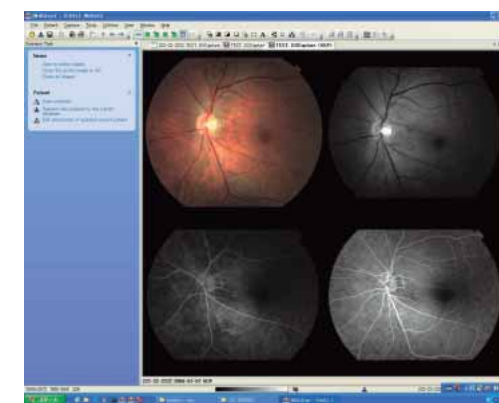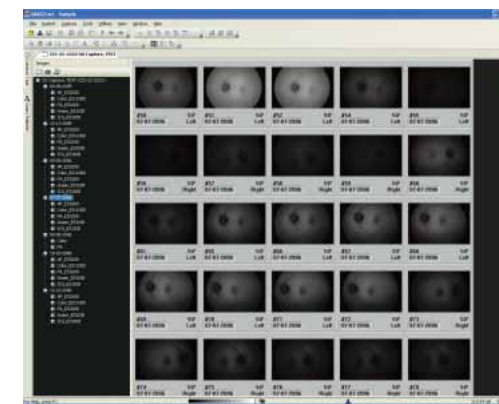


PDT

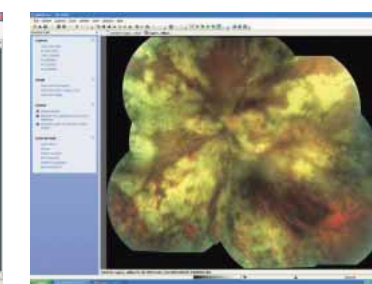

Auto Mosaic



Patient Consultation

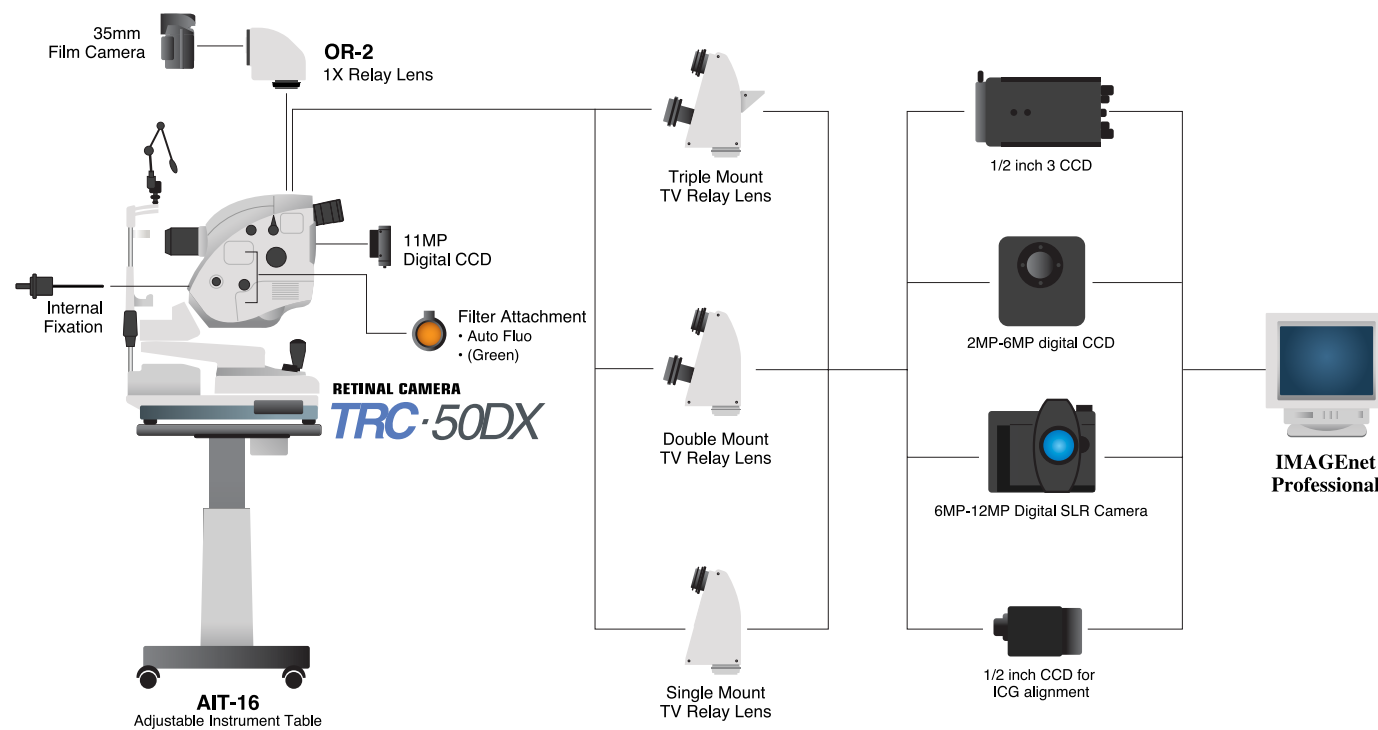

*All cameras are other manufacturer's products.

"In compliance with the terms of the Export Administration Regulation of the United States of America, some optional software modules may not be available in all regions or countries."

Specifications

	TRC·50DX	TRC·50DX (Type 1A)
Capture Mode	Color/Red-Free/FA	Color/Red-Free/FA/ICG/Auto Fluo
Auto Fluorescence	OPTION*	○
Angle of Coverage	50°/35°/20°	
Photographic Magnifications	With 35mm camera at zero diopter ;1.84x at 50°, 2.45x at 35°, 4.28x at 50°	
Working Distance	39 mm	
Diopter Compensation Range for patient's eye	0 -10D~+6D	
	- -23D~-9D	
	+ +5D~+23D	
	A +22D~+41D	
Diopter Compensation at finder	-6D~+5D	
Split Focus	○*	
Internal Fixation	OPTION*	○
Counter	○*	
Optical Head tilt	UP 15° / DOWN 10°	
Light Source	For photography : Max 300WS Xenon	
Power supply	100-120 V,200-240V / 50-60 Hz	
Power Consumption	1500VA	
Dimensions / Weight	340 mm(W) x 505 mm(D) x 506-715mm(H) / Aprox. 35 Kg	

* is available in not all countries



*Subject to change in design and/or specifications without advanced notice.

IMPORTANT

In order to obtain the best results with this instrument, please be sure to review all user instructions prior to operation.



TOPCON CORPORATION



TOPCON CORPORATION



EC9711081



TOPCON CORPORATION

75-1 Hasunuma-cho, Itabashi-ku, Tokyo 174-8580, Japan. Phone:3-3558-2523/2522 Fax:3-3960-4214 www.topcon.co.jp

TOPCON MEDICAL SYSTEMS, INC.
37 West Century Road, Paramus, NJ 07652, U.S.A.
Phone: 201-261-9450 Fax: 201-387-2710 www.topcon.com

TOPCON CANADA INC.
110 Provencher Avenue, Boisbriand QC J7G 1N1, CANADA
Phone: 010-45865077 Fax: 010-2844940 www.topcon.eu

TOPCON EUROPE MEDICAL B.V.
Essebaan 11, 2909 LJ Capelle a/d IJssel, THE NETHERLANDS
Phone: 010-45865077 Fax: 010-2844940 www.topcon.eu

ITALY OFFICE
Via De' Industria n. 60, 20037 Paderno Dugnano, (Milano), ITALY
Phone: 02-9186671 Fax: 02-91081091 www.topcon.it

DENMARK OFFICE
Praestemarksvej 25, 4000 Roskilde, DENMARK
Phone: 046-32-7500 Fax: 046-32-7555

IRELAND OFFICE
Unit 69 Western Parkway Business Centre,
Lower Ballymount Road, Dublin 12, IRELAND
Phone: 01460-0021 Fax: 01460-0129

TOPCON S.A.R.L.
89, Rue de Paris 92585 Cligny, Cedex, FRANCE
Phone: 01-41069494 Fax: 01-47390251 www.topcon.fr

TOPCON DEUTSCHLAND G.m.b.H.
Giesserallee 31-33, 47877 Willich, Germany
Phone: 02154-8850 Fax: 02154-885111 www.topcon.de

TOPCON SCANDINAVIA A.B.
Neogatan 2, 43121 Mölndal, SWEDEN
Phone: 31-7109200 Fax: 31-7109249 www.topcon.se

TOPCON ESPAÑA S.A.
Frederic Mompou 5, Edificio Euro 3, 08960 Sant Just Desvern Barcelona, SPAIN
Phone: 34-93-4734057 Fax: 34-93-4733932 www.topcon.es

TOPCON (GREAT BRITAIN) LTD.
Topcon House, Kermet Side, Bone Lane, Newbury, RG14 2PX U.K.
Phone: 01635-551120 Fax: 01635-551170 www.topcon.co.uk

TOPCON SOUTH ASIA PTE.LTD.
Blk 192 Pandan Loop, #07-01 Pantech Industrial Complex, Singapore 128381
Phone: 62780222 Fax: 62733540 www.topcon.com.sg

TOPCON INSTRUMENTS (MALAYSIA) SDN. BHD.
Exeña Business Park Block C, 1st Floor, Jalan Ampang Putra,
Taman Ampang Hill, 55100 Kuala Lumpur, MALAYSIA
Phone: 03-42701192 Fax: 03-42704508

TOPCON INSTRUMENTS (THAILAND) CO.,LTD.
77/162 Simsathorn Tower, 37th Floor, Krungthoburi Road,
Klongsai, Klongsam, Bangkok 10600 THAILAND
Phone: 02-440-1152-7 Fax: 02-440-1158

TOPCON KOREA CORPORATION
2F Yooeung Bldg., 159-3, Seocho-Dong, Seocho-Gu, Seoul, 137-876, Korea
Phone: 2-2055-0321-7 Fax: 2-2055-0319 www.topcon.co.kr

TOPCON OPTICAL(H.K.) LTD.
2/ F., Meeco Industrial Bldg., No.53-55 Au Pui Wan Street,
Fo Tan Road, Shatin, N.T., Hong Kong
Phone: 852-2690-1328 Fax: 852-2690-2221 www.topcon.com.hk

TOPCON CORPORATION BEIJING OFFICE
Building B No.9, Kangding Street
Beijing Economic Technological Development Area, Beijing, 100176 CHINA
Phone: 10-6780-2799 Fax: 10-6780-2790

TOPCON CORPORATION SHANGHAI OFFICE
14L HUAMIN EMPIRE PLAZA, No.726, Yan-an Xi Road, Shanghai, CHINA
Phone: 21-52387722 Fax: 21-52370761

TOPCON CORPORATION BEIRUT OFFICE
P.O. BOX 70-1002 Antelias, BEIRUT-LEBANON
Phone: 961-4-523525/523526 Fax: 961-4-521119

TOPCON CORPORATION DUBAI OFFICE
C/O Atlas Medical FZCO., P.O.BOX 54304, C-25,
Dubai Airport Free Zone, United Arab Emirates
Phone: 971-4-2995900 Fax: 971-4-2995901